


POWERSHAPE 2™ RF / Vacuum / Low level laser

The RF device for professionals is used for body contouring and face lifting. The handpiece of Powershape2 combines three technologies such as RF, vacuum, and low-level laser, and their simultaneous operation is expected to bring about synergetic effect in terms of obesity improvement. In addition, roll handpiece mounted with motorized dual rollers as well as multi polar handpieces of various sizes are provided to allow a sophisticated care.

Structure of handpiece


Roll handpiece

Vacuum-assisted motorized roll handpiece enables both Laser & Bi-polar RF energy precisely targeted where it can have maximum effect.

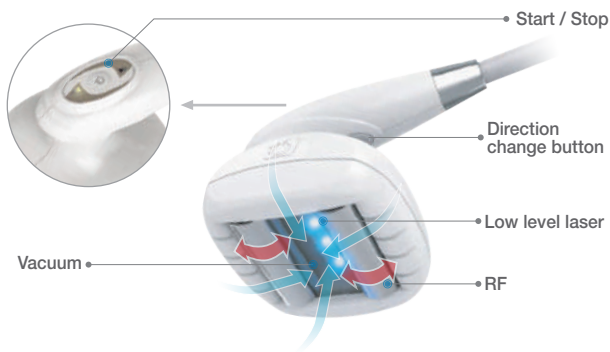


Multi-polar handpiece

Multi-polar RF technology painlessly delivers focused and optimal RF energy to the skin surface and to subcutaneous fat layers.

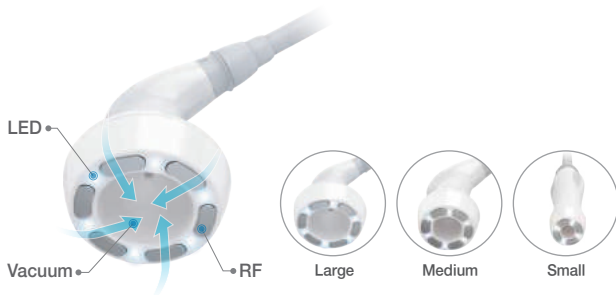

Thermal handpiece

Increased thermal energy conductivity that reduces treatment time thanks to the ceramic coated bi-polar tips.

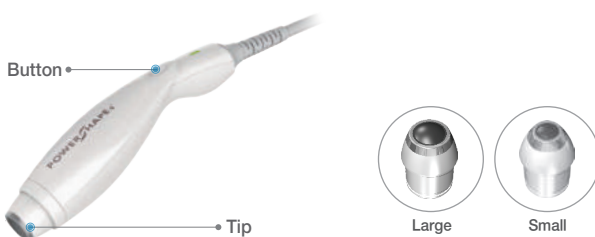

Specification

Energy type	RF / LLL / Vacuum
Frequency	2MHz
Vacuum pressure	300mmHg
Wavelength	660 nm
Input Power	AC 100~240V, 50 / 60 Hz
Dimension (WxDxH)	414 x 523 x 1,230 (mm)